



Q5k

DIALYSIS BED

Registration code: SKD-C




SAIKANG / Medical

Healthier
your life

Q5k

FEATURES AND HIGHLIGHTS

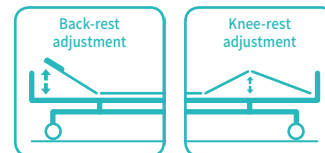

Technical parameters

▪ External size (LxWxH):	1930x880x600 mm
▪ Mattress platform:	1710x780 mm
▪ Safe working load:	200 kg
▪ Back-rest adjustment:	0-75° (±5°)
▪ Knee-rest adjustment:	0-45° (±5°)

Technical configuration

<input checked="" type="checkbox"/> Linear motor	2pcs
<input checked="" type="checkbox"/> Aluminium alloy foldable side rail	1set
<input checked="" type="checkbox"/> PP bed ends with safe lock	1set
<input checked="" type="checkbox"/> PP bed platform, removeable, washable	5pcs
<input checked="" type="checkbox"/> 5" double side castors	4pcs
<input checked="" type="checkbox"/> IV pole	1set
<input checked="" type="checkbox"/> IV pole prevision	4pcs
<input checked="" type="checkbox"/> Drainage hook	2pcs

QUALITY CONTROL

SAIKANG medical, as a pioneer in the field of medical furniture, always choose branded suppliers strictly to ensure the safety and reliability of raw material. SAIKANG adopt robot welding, automatic epoxy coating, efficient assembly line, dedicated operator, rigorous quality control, only to guarantee SAIKANG's commitment on quality



Design

This Bed uses unique design to meet hospital strict gremfree request. Head & Foot board and platform are all seamless, which supports effective sterilization.



Steels | welding

Baosteel from Fortune Global 500. Panasonic Robotic ensure 360° full smooth welding.



Plastic | Painting

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.
11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.



Electric Motor

Linear motor,with UL /ROHS/EN standard.Noise under 50 dB,with lifetime 20,000 times movement.

Accessories



Platform

Five sections PP Platform with Integrated mattress Retainer, removable for sufficient sterilization.



Castors | Brake padel

5" Double covered castors with TPR tire, no worn-out after running 30KM, united forming without bolts

Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time.

The double-folded large-area rubber design, horizontal and vertical two-sided rubber, allows the nurse to touch the brakes when the soles of the feet are soft and large-area plastic surface, saving effort and protecting the upper.


Head & Foot Board

Modern design PP head&foot board.
Easy to be cleaned and installed.



Handheld controller

Handheld controller for patients and nurses, ergonomic design.

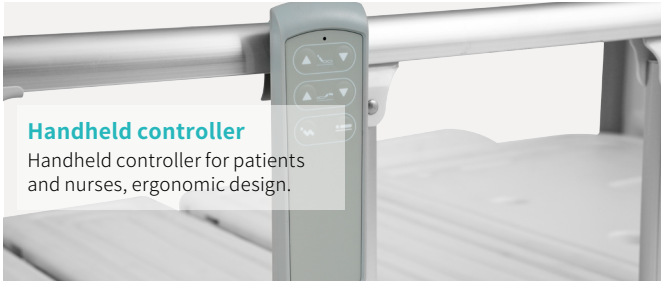


Siderail

Side rail made from aluminium alloy and take 4-staff design. And it with aluminium alloy protection for patients can get on or off bed easily.



Drainage Hooks

Mobile plastic drainage hooks.

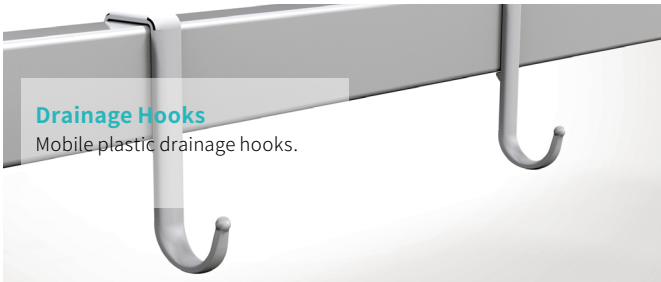

IV Pole

Full SS IV pole, bearing 2.5kg per hook.



Safe working load
200kg



Back-rest adjustment
75°



Knee-rest adjustment
45°



CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service

WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc
 W3400 Corrientes, Argentina
 TEL: +5493794312231
 Av. Senador Vergueiro, 4204 - sala 31
 Rudge Ramos - São Bernardo do Campo - SP, Brazil
 TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd
 Mezzanine Floor, Taj Tower,
 upper hill, nairobi, Kenya
 TEL: +254714084806

South East Asia Branch

Saikang Medical Technology Ltd
 Sala Thammassop Subdistrict
 Thawi
 Watthana District
 Bangkok Thailand
 TEL: +66 642741046

UK Branch

Uk Healthward International Limited
 DEPT 111,
 196 High Road Wood Green,
 London, N22 8HH
 United Kingdom
 TEL: +31 (0) 202111106

Russia Branch

LLC 'SAIKANG'
 454004, Russian Federation,
 Chelyabinsk,
 Universitetskaya Naberezhnaya,
 Russia
 TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC
 Musaffah
 P.O Box 53200, Abu Dhabi,
 UAE
 Tel: +971 585710652

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
 STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue